

 **powerhands**  
by **POWERFIT**



# User Manual

Read all instructions before using and keep this manual for future consultation.



# Introduction

Thank you for using Powerhands, enjoy a pleasant hand massage anytime and anywhere.

In order to use and maintain this product correctly, please read the manual carefully before using, and keep it for further inquiry.

## Table of content

Introduction .....	3
Safety Information .....	4
Charging .....	4
Control panel .....	5
Control panel functions .....	5
Troubleshooting .....	7
Technical information .....	7

## Safety Information

Please read this instruction carefully before use.

This product is suitable for people who use a computer for a long time, play with mobile phones, suffer from finger numbness and want a relaxing massage on their hands.

This product needs to be charged by USB with a mobile phone charger. Please note that only chargers with output voltage 5V 1Amp - 5V 4Amp can be used, and no chargers with output voltage greater than 6V can be used. The red color indicator will flash during charging and the indicator will turn green after fully charged.

Turn off the device to depressurize the interior before removing your hand.

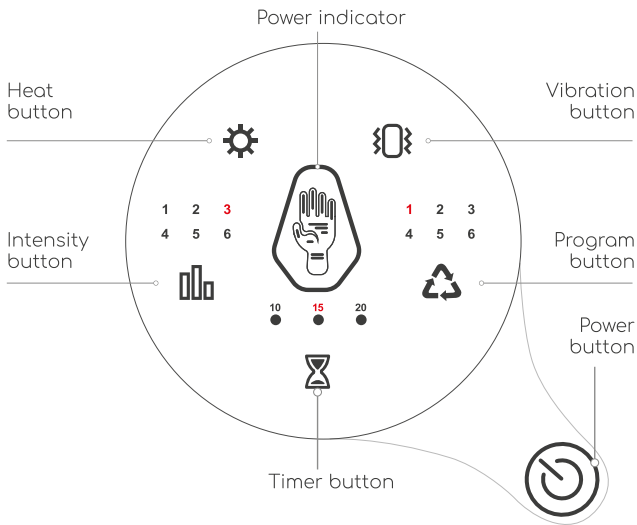
## Charging

If the Powerhands battery is depleted, it will need near to 4 hours to fully charge it. Use the USB cable provided with a wall charger (not included).



Note that Powerhands does not support charging by computer.

Once is fully charged, with vibration and heating on, Powerhands can work no less than 100 minutes.

## Control panel



## Control Panel Functions

-  **Power button:** Turn Powerhands on or off. Press for 2 seconds to turn it on, it will start to work with the default setting mode: P1 mode, level 3 intensity and timer in 15 minutes. To turn it off press 2 seconds again.
-  **Heat button:** Press it to turn on/turn off the heating function. The symbol changes to red, it would take up to 5 minutes to reach the maximum temperature.



## Battery level indicator

100 - 51%



50 - 21%



20 - 1%



When the estimated remaining battery charge drops down to approximately 1%, the device will beep for 4 seconds, and then it'll shutdown in 1 minute.



**Intensity button:** The default value is L3, with a total of 6 levels of intensity. You can select it from L3-L4-L5-L6-L1-L2 in sequence L1 is weakest and L6 is strongest.



**Timer button:** The default startup time is 15 minutes. Press the Hourglass symbol to switch between 10 minutes, 15 minutes or 20 minutes, it can be adjusted freely.



**Vibration button:** Press to turn vibration on/off. There are 2 vibration modes:

- **Pulse type** (symbol changes to blue color)
- **Continuous vibration** (symbol changes to red color), you can select the vibration by touching the symbol.



**Program button:** The default program is P1, with a total of 6 programs available. You can choose your favorite program to massage, which is from P1-P2-P3-P4-P5-P6-P1 in sequence.

### Programs description:

**P1:** Combination of P2 to P6. The lead time of each mode is 20 seconds.

**P2:** Palm + Wrist randomly

**P3:** Palm + Wrist simultaneously

**P4:** Wrist

**P5:** Palm + Wrist alternately

**P6:** Palm

Powerhands comes with memory function, It will record the setting of last use.

## Troubleshooting

Fault Type	Fault Cause	Fault solution
Unable to start product	The battery has run out of power	Recharged
	Intensity mode is not enabled	Turn on Intensity mode
Make a faint sound when used	Sound of the air pump (Buzzing)	Machine running sound, no trouble
	Solenoid valve sound (tick-tock)	

## Technical information

Product Name	Powerhands
Model name	R-H001
Voltage	DC 5V
Power	5W
Charging	USB Type-C Interface 5V
Battery Capacity	3.6V, 2500mAh Lithium Battery
Suggested time	15 minutes
Heating Function	Up to 43° C
Weight	990g.

## Warranty

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.



## Included items:

- 1 Powerhands massager
- 1 User manual
- 1 USB-C cable
- 1 Travel bag
- 1 Finger massager



Made in China. ©2022. All rights reserved.

### DISTRIBUTED BY:

Orange Planet  
De Els 3a  
5107 RK Dongen  
The Netherlands



Download the manual  
in different languages here

[www.yellowweb.nl/nl/manuals-powerhands](http://www.yellowweb.nl/nl/manuals-powerhands)



**LE TRI**  
**+ FACILE**

**EMBALLAGES A SEPARER**  
**ET A DEPOSER DANS LE**  
**BAC DE TRI**



*Séparez les éléments avant de trier*